ELECTRO BF

PL

INSTRUKCJA

EN

USER MANUAL

CZ

NÁVOD

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







Programy

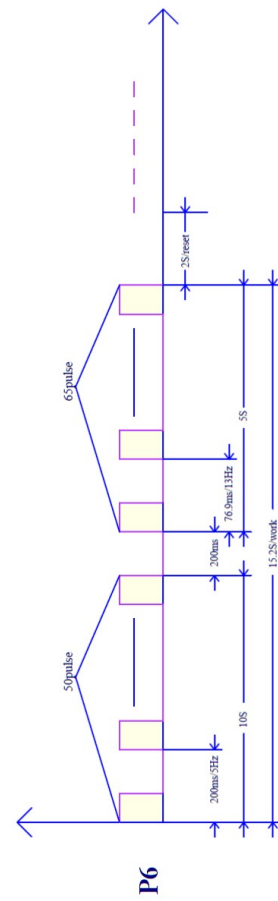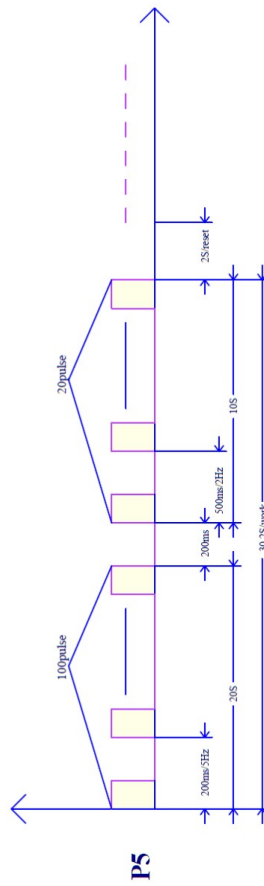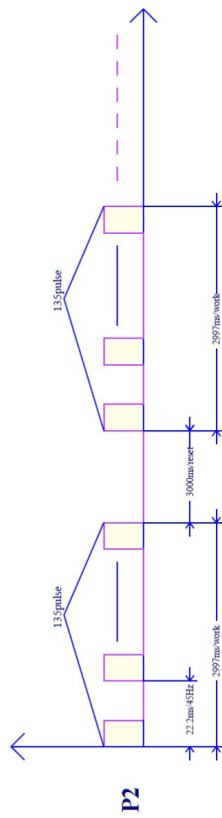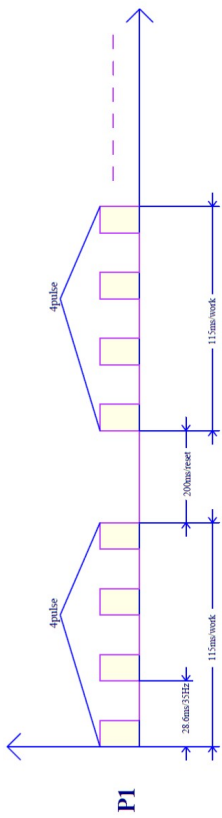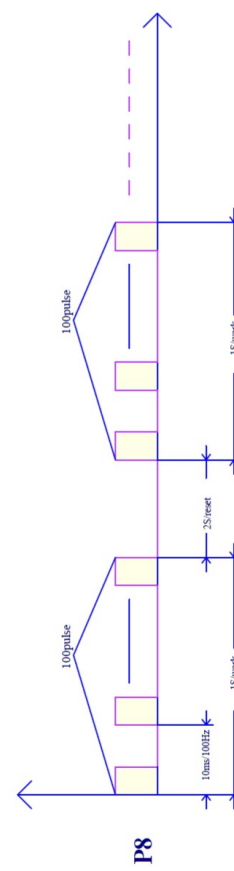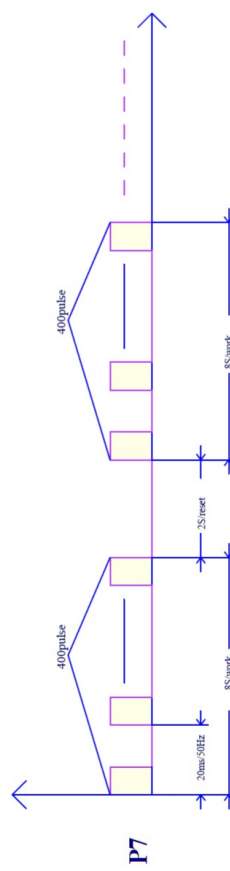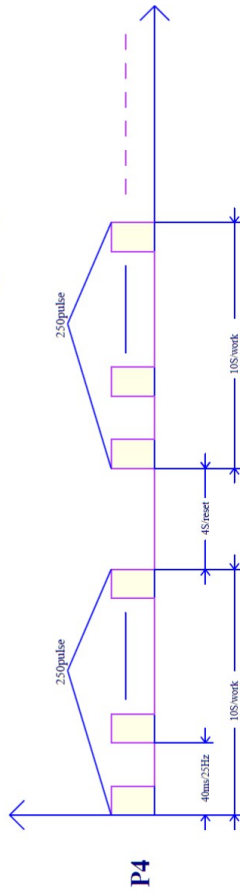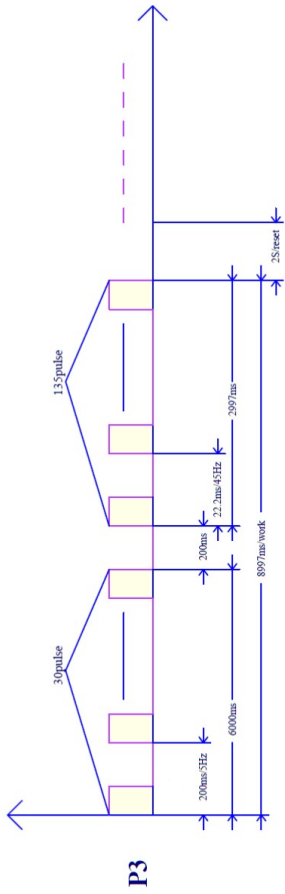**EN**


Precautions: Study the precautions and instructions included in this manual before using the equipment. Safety instructions are provided for your own safety in order to avoid harm to other persons as well as damages to property. Following these recommendations ensures correct and safe use of the product.

Warning:

1. **Do not use the massager if you use the following medical devices because it can cause injuries of your body and damage to the equipment.**
 - Life support devices, i.e. pacemaker and similar devices
 - Medical equipment such as heart valves and similar devices
 - Electrocardiogram (ECG) and similar devices
2. **It is forbidden to use the equipment in the following cases, since its use can lead to injuries and skin damage**
 - Patients with heart diseases
 - Pregnant women and postpartum women
3. **Consult your physician before using the massager if any of the following health problems applies to your situation:**
 - Severe allergies
 - Patients with temperature over 38°C
 - Skin problems
 - People during recovery
 - People suffering from acute or long-lasting diseases, bone fractures, sprains, etc.
 - People having problems with blood clotting
 - People with osteoporosis and cancer
 - Epileptics
 - People with a contraindication to use hydrogel
 - Diabetics
 - People in the course of medical treatment monitored by doctors, etc.
4. **Do not use the equipment on the following body parts since it can cause injuries and other damages**
 - Chest, arms, neck and their vicinity (head, face, etc.)
 - Injured body parts, intimate parts
 - Curing scars, on medical devices placed under your skin
 - On tattooed skin
5. **Proceed according to the following instructions in order to avoid damages, injuries, etc.:**
 - Do not put the device into a microwave, oven, do not bring it close to fire
 - Do not use the equipment against its intended purpose
 - Do not use it on animals and children
 - Do not put any lubricates and liquids on the device
 - Do not use it when driving or exercising
 - Do not use the device in the bathroom and rooms of increased temperature
6. **The product is intended for adults only. The product should be kept beyond the reach of children. Contact a doctor immediately if the battery is swallowed by a child.**
 - Do not dismantle, repair and modify the device on your own – you may damage the equipment as well as bodily harm.
 - Protect the equipment against falling or damage.
 - The use of damaged equipment can result in injuries.
 - Avoid the contact of electric transmitters with metal parts – it can lead to electric shock.
 - Do not use the massager when using other electronic devices since it can result in undesired side effects.
 - Stop using the product immediately and consult your physician if any of the following symptoms appear when you use the product: pain, discomfort, skin dryness, convulsions, burning sensation or other worrying symptoms. Ignored problems can get more serious.
 - Switch off the device and remove the batteries immediately after you notice unpleasant smell, water inside electronic parts or any disturbing sounds.
 - For hygienic reasons the device shall be used by one person only
7. **Check the condition of batteries included in the sets on regular basis – damaged batteries can cause injuries or damage the equipment**
 - Switch off the device and remove the batteries immediately after you notice unpleasant smell, water inside electronic parts or any disturbing sounds.
 - For hygienic reasons the device shall be used by one person only
 - If the liquid leaks from the battery, remove the battery immediately from the device and contact the distributor
 - If the liquid from the batteries gets into your eyes, immediately wash your eyes with a lot of water. And contact your physician. Rubbing your eyes can lead to their injury.
 - If the liquid from the batteries comes into contact with your skin, immediately wash your skin with a lot of water.
 - Make sure that the battery is installed according to symbols on the device
8. **Check the condition of hydrogen-covered pads on regular basis – damaged pads can cause injuries**
 - Do not use damaged or dirty pads
 - Do not bend the pads in order to avoid their damage
 - Make sure that the pads stick closely to your skin. If the hydrogel does not touch the skin with its all surface, electrical impulse can be concentrated in one spot only.
 - Switch off the device before taking it off.
 - Make sure that the whole surface of the electrode is covered with gel. Damaged spot (where the electrode is not covered with gel) should not contact your skin

List of parts from Electro BF:

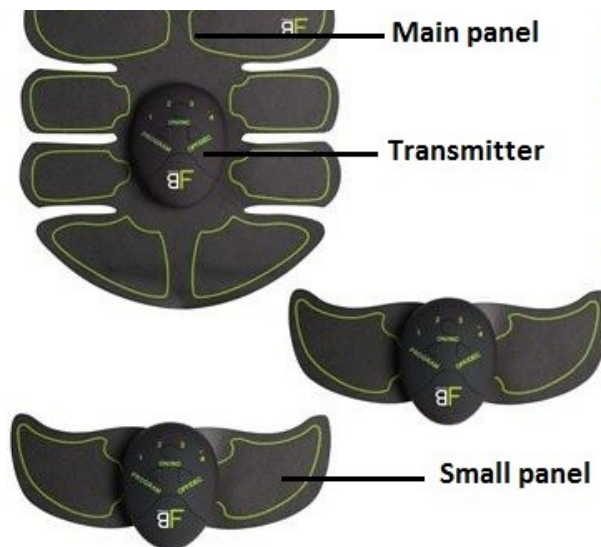Set includes:
Panel: 3 pcs.
Transmitter: 3 pcs.
Use manual: 1 pc.

Before using the equipment, put batteries (2*AAA) in transmitter

Note: Batteries are not included

Keys:

Warning:

- The device should be switched on only after it has been correctly placed on the body and it should be switched off before it is taken off. After you finish using the product, detach the transmitters from stimulating pad
- Immediately stop using the product and consult your physician if any of the following symptoms appear when you use the product: pain, discomfort, skin dryness, convulsions, burning sensation or other worrying symptoms.
- Excessive and long-lasting use of the device can lead to reddening of skin and exhaustion of muscles
- Take off jewellery before starting your exercise
- Remove your watch and put other electronic equipment at a safe distance from your body in order to avoid electric shock and short circuit in the devices
- Do not use the device on skin covered with cream or other lubricants – it can reduce grip of the pads with conductive gel
- The pads will stick better to hair-free spots – however, the device cannot be used on irritated skin on the day the hair has been removed
- Do not heat and cool the equipment – it can result in its deformation and damage
- Before using the massager, wash carefully the part of your body where it will be used
- The product is not waterproof.

Batteries:

Proceed according to the following instructions concerning batteries

- If the device is not used for a long period of time, remove the batteries.
- Do not use outdated batteries.
- Do not touch the batteries with wet hands.
- Use AAA, CR2032 batteries or battery packs
- Batteries wear out faster in rooms with low temperature
- The longer the device is used, the weaker impulses will be. This is caused by reduced battery capacity. If you notice that the device is not working with the same output as when you started using it, replace the batteries or charge them.

Battery replacement:

1. **Open the protective cover on the back side of the transmitter**
2. **Insert the battery in its proper place with “+” side up**
3. **Insert the protective cover**

How to use conductive gel:

- Do not use gel pads that have visible damages, are worn-out or dirty,
 - The same procedure applies to small and big pads
 - If the pads are incorrectly attached, the device will not work properly
1. **Conductive gel wears out during the product use. Replace a gel pad when:**
 - The gel loses its viscosity
 - The gel quality is significantly worsened
 2. **When using the device:**
 - store the pads in a dry and cool place, away from sunlight and high temperatures which might destroy the gel
 - do not put pads on dirty surface. Dirty pads can cause discomfort during exercise
 - Depending on the frequency of use, skin condition and storage, the gel can lose its viscosity after it is used 25-30 times. In such case, follow the procedure for cleaning the device.
 - If the gel sticks to your clothes, do not tear it off. Wet the fabric slightly and gently pull the gel.
 3. **Replacement of a gel pad**
 - I. Remove dirty or damaged pad
 - II. Install a new pad

Method of use:

1. **Take the product out of its case (if included)**
2. **Attach a transmitter**
3. **Remove protective foil from the device**
4. **Place the device on the body**

Intensity

The device offers 6 work modes and 10 levels of intensity.

NOTE: All 3 pads can be used at the same time.

Device removal

- Make sure that the device is switched off and then pull it gently to take it off.
- Try not to touch gel with your fingers.
- Put protective foil on the device and put the device in its case. If the gel pads are dirty, replace them.

1. After completed your training, gently detach the pads from your body
2. Put protective foil on the device and put the device in its case
3. Switch off the device before taking it off.

If there are gel traces on your skin, wipe them off with water and soap. If the gel remains on your skin too long, the skin might get irritated.

Correct way to take off the device: Grab the top part of the device and detach it gently

Incorrect way to take off the device: Grabbing one corner and pulling the pad by a part

Maintenance and storage

1. Place the device on protective foil after completed workout.
 - Do not fold gel pads or glue them together
 - Do not press hard the gel pads to protective foil
 - Do not fold the protective foil
2. Put the device into its case after completed workout
 - Make sure the protective foil and the case are clean
 - Avoid exposure to high temperature
 - Use the device only according to its intended purpose
 - Store the device far away from fire
 - Remove the batteries from the device if it is not used for a longer period of time (about 1 month)

Disposal

The device shall be disposed according to local laws on the territory where it remains. In case of long-time use, fast wearing or damaged parts should be replaced and disposed.

Make sure that the device is switched off before its maintenance.

- Do not wash the device with water.
- Do not use removers to clean the device
- Wash your hands before touching gel pads
- Use a slightly wet cloth for cleaning in order to avoid damaging or scratching the device.
- If the device is very dirty, use a water-wetted cloth with gentle washing liquid to clean it.
- After cleaning, put the device in a cool, dry area.
- If gel is dirty, it will be difficult to stick it to your body. In such case, proceed according to the following steps:
 - Put a few drops of water on your finger and smear them on the gel pad to remove dirt
 - Next, leave the pad to dry in a cool and dry area. Put it with the gel side up
 - Put protective foil on the device when the gel is dry









Possible problems and troubleshooting

Problem	Reason	Repair procedure
The device does not switch off	There are no batteries in the device	Insert the battery
	Battery is not correctly installed	Put the battery with "+" up
	Battery is empty	Replace the battery
	Battery is worn out	
	Incorrect battery	Use battery CR2032
I don't feel any impulse	Gel does not stick to your skin	Make sure there is no gap between the gel and your skin
	Gel pad is damaged	Replace the gel pad
	Gel is dirty	Clean the gel pad
Impulse is too weak	Intensity level is too low	Increase the impulse intensity
	Low power level	Replace the battery
	Gel pad is damaged	Replace the gel pad
	Gel is dirty	Clean the gel pad
	Intensity level is too high	Reduce the intensity level
Impulse is too strong	Gel does not stick to your skin	Make sure there is no gap between the gel and your skin
	Battery is empty	Replace the battery
The device stops working suddenly	Battery is empty	Replace the battery
Problem	Reason	Repair procedure
The device does not stick to the body	Gel pad is dirty	Clean the gel pad
	Gel pad is damaged	Replace the gel pad
	Gel pad is not glued correctly to the skin	Make sure there is no gap between the gel and your skin
Skin turns red and irritated	Gel pad is dry and dirty	Clean the gel pad
	Gel pad is damaged	Replace the gel pad

If the above suggestions are not applicable, make sure that there is a required point of contact with the battery and press the protective flap. Gel is gets worn out during exercise – buy a new pad if needed. If the problems with the device continue, stop using it and contact distributor.

Technical specification of Electro BF:

Weight	57.4 kg
Dimensions	210x100x23 mm
Output	Max. 9.8mA
Frequency	1~ 100 Hz
Voltage	Alkaline 1.5V AAAX2EA

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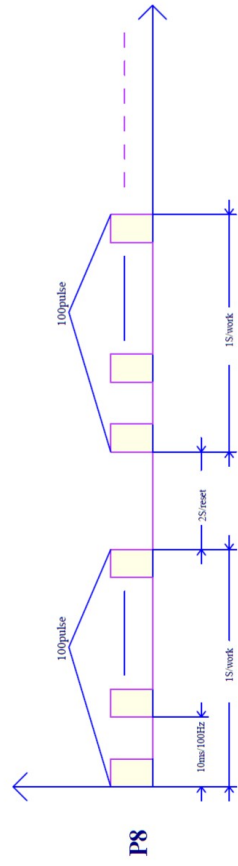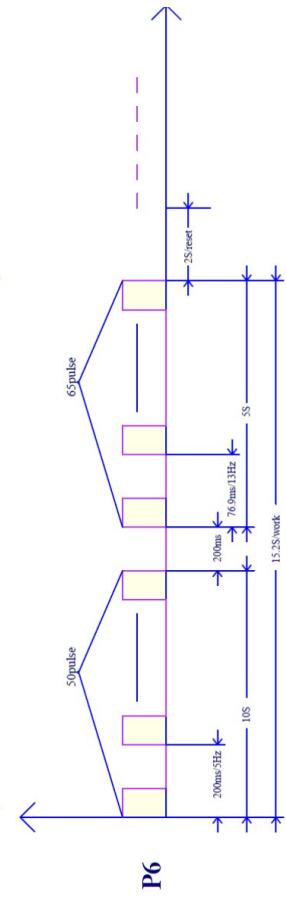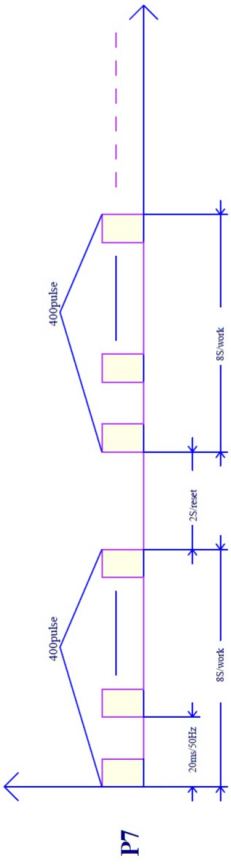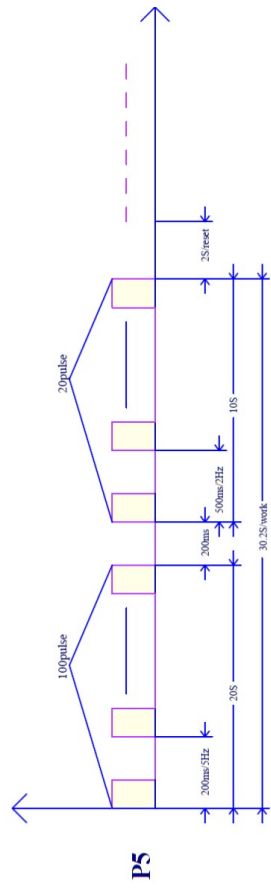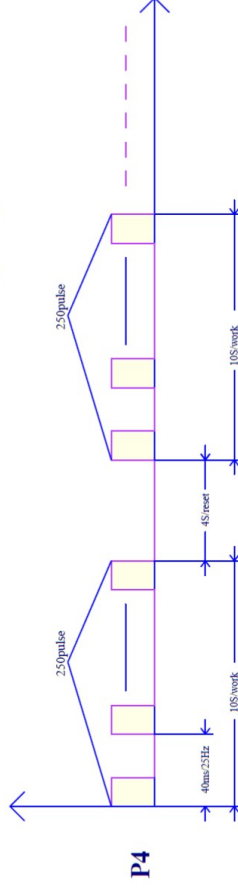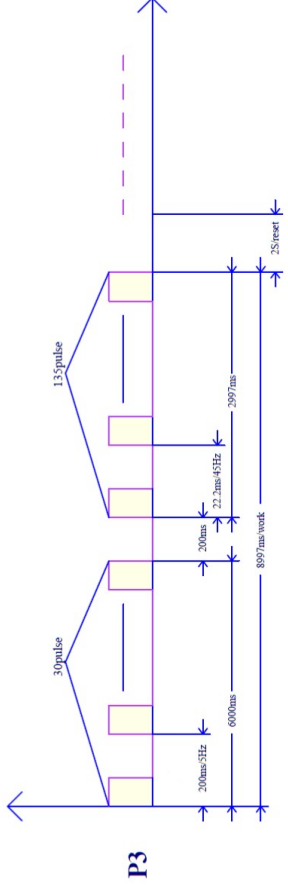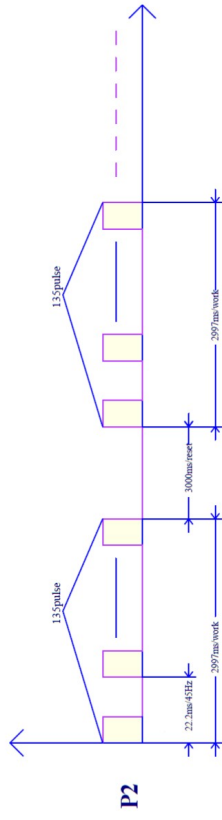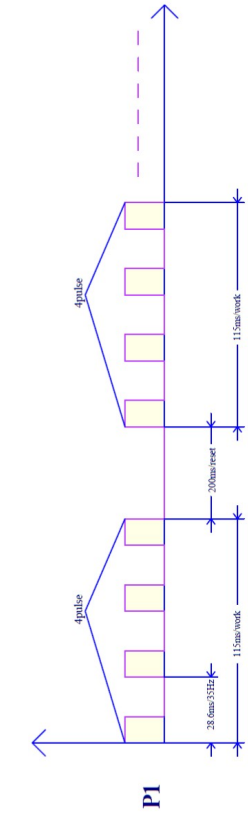EN: The marking of the equipment with the crossed-out waste container symbol indicates that it is prohibited to place used electrical and electronic equipment with other waste. According to the WEEE Directive on the management of used electrical and electronic waste, separate disposal methods must be used for this type of equipment.

The user who intends to dispose of this product is obliged to take it to a collection point for used electrical and electronic equipment, thus contributing to reuse, recycling or recovery and thus to the protection of the environment. To do this, please contact the point from which the appliance was purchased or your local authority. Hazardous components contained in electronic equipment may cause long-term adverse effects on the environment, as well as harmful effects on human health.

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Date of sale:

GUARANTEE TERMS:

1. The Seller, on behalf of the Guarantor, provides a guarantee in the territory of the Republic of Poland for a period of 24 months from the date of sale.
2. The guarantee will be honoured by the shop or service centre upon presentation by the customer:
 - a legibly and correctly completed warranty card with the sales stamp and the seller's signature,
 - a valid proof of purchase of the equipment with the date of sale / receipt, the goods claimed.
3. Any defects and damage discovered during the warranty period will be repaired free of charge within a maximum of 21 days from the date of delivery of the goods to the service.
4. In the case of the necessity to import parts, the repair period may be extended by the time necessary for their import, but not longer than 90 days.
5. The warranty does not cover:
 - mechanical damage and defects caused by them,
 - damages and defects resulting from improper use and storage,
 - improper assembly and maintenance,
 - damage and wear of components such as cables, straps, rubber parts, pedals, sponge grips, wheels, bearings, etc.
6. The warranty is void in the event of:
 - expiry date,
 - self-repair,
 - failure to observe the rules of correct operation.
7. Product returned for repair should be complete and clean. In the case of defects the service has the right to refuse acceptance for repair. If the product is delivered dirty, the service centre may refuse to accept it or clean it at the customer's expense with his written consent.
8. The warranty does not cover installation and maintenance work, which, according to the user manual, must be carried out by the user himself.
9. The guarantor also informs that it provides post-warranty service.
10. The goods should be protected for shipping.
11. In order to make use of the warranty, please follow the procedure on the website: <https://serwis.abisal.pl/>.

In case of non-conformity of the sold thing with the contract, the buyer is entitled by law to legal remedies from and at the expense of the seller. The guarantee does not affect such remedies.

THE EQUIPMENT IS NOT INTENDED FOR REHABILITATION AND THERAPY

NOTES ON THE COURSE OF REPAIRS

Item	Date of notification	Date of provision	Course of repairs	Signature of the recipient (shop, owner)

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